


# Guide for securing an oral endotracheal tube



**Assemble equipment**




**Apply Cavilon around lips and cheeks**

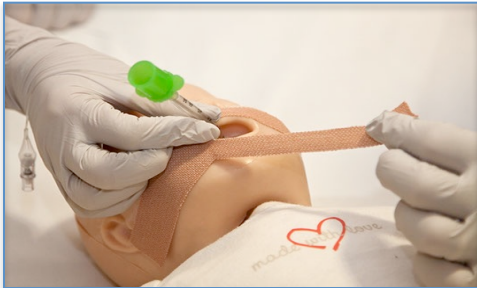

**Place V of trouser leg at angle of mouth. The upper part of the trouser leg is then fixed along the upper lip**




**The lower part of the trouser leg is then turned back to wind around the ETT in an upward spiral**

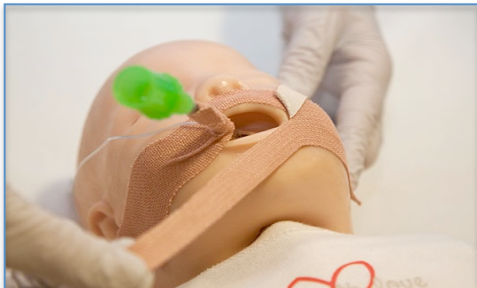


**Using the second trouser leg, place the V at the opposite angle of the mouth. The lower part of the trouser leg is then fixed below the lower lip**



**The upper part of the trouser leg is then fixed above the upper lip**

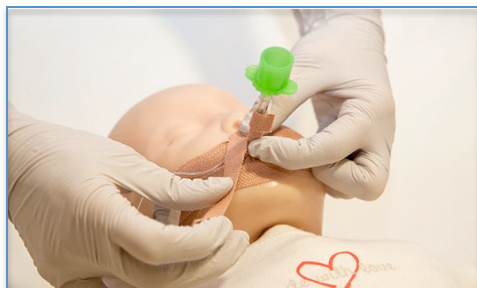


**The free end of the trouser leg is then wound over the ETT in a downward spiral**

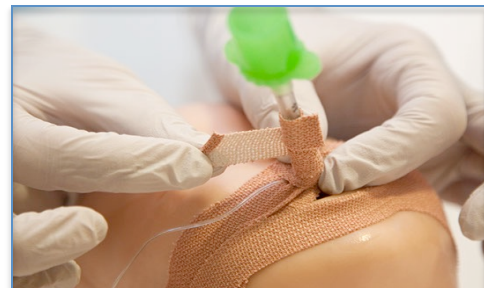

**Fold the last part of ET tape back on itself to make a small flap to facilitate removing the tape later**



**The third piece of tape has a central slit through which the ETT connector is passed**

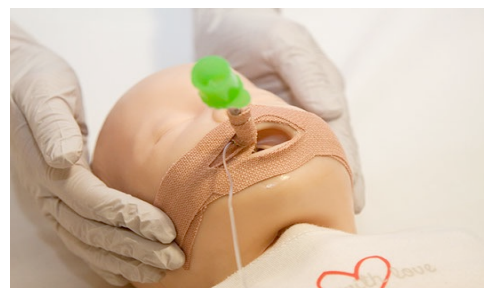

**This is then applied above and below the lips**